

# Filter Slider Pieces

quantity: one of each


Steve Smith

Bustamante Group

phone 666-2761


steve@alice.berkeley.edu

outer slide

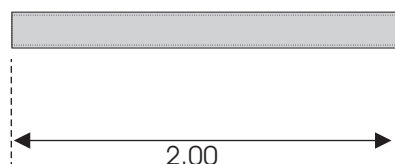



square telescoping brass tubing  
3/16" outside, 5/32" inside

inner slide

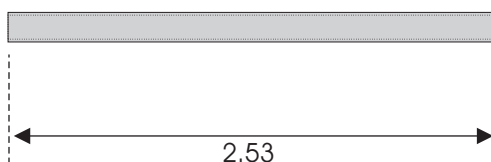



square telescoping brass tubing  
5/32" outside, 1/8" inside

inner cable guide

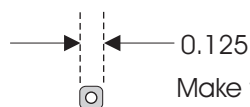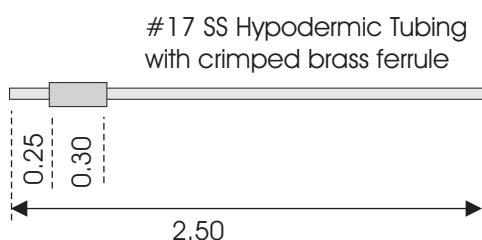



Make ferrule 0.140 round and drill #53.  
Then place on hypo tube and whack sides  
with hammer and file down until it  
is square and fits snug inside  
inner slide tube (1/8" x 1/8")  
Ferrule will be soft-soldered into inner slide.

handle cable guide

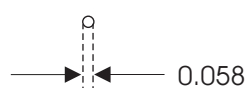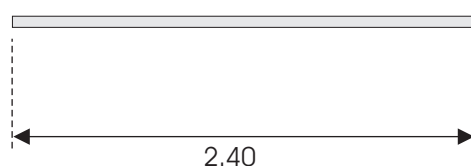

#17 SS Hypodermic Tubing

slider handle

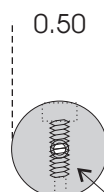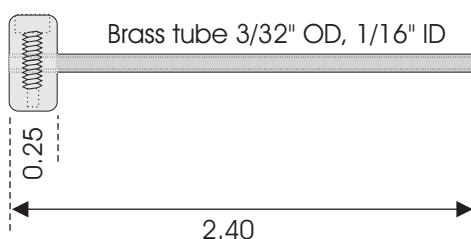

Press-fit brass handle tapped and  
counterbored for 4-40 x 3/8" Allen cap

1/32" steel cable  
comes out here

screw clamps end of  
cable inside knob